




# TELE-VIEW® DYNAMIC EQUINE EXERCISE ENDOSCOPE

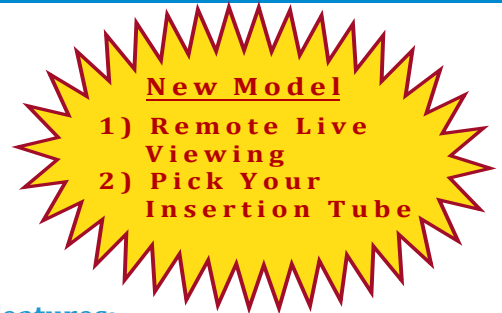

Control Body Gives 90° Up/Down Articulation!

Irrigation Bulb

Android Phone Records Audio & Video

TV-507/508 Model

### Features:

- » High Resolution Images and Video
- » Records Audio and Video to Android Phone
- » Live Viewing Capability
  - Easily Mirror Android Image up to 10 meters to Another Phone or Tablet
  - Watch Remotely with Strong Cellular Signal on Phone, Tablet or Computer

### Benefits:

- » Simple and Easy Design
- » 'Gold Standard' for Diagnosing Airway Restriction Problems
- » Results in 30 Minutes
- » Well Tolerated by the Horse

### Available Insertion Tubes

Model TV-507  
11.5mm Stabilized



Model TV-508  
9.8mm Flexible

The Tele-View® Dynamic Equine Exercise Endoscope is a 'simple and easy' to use method for diagnosing upper airway breathing problems. Dynamic endoscopy has been proven to be the 'Gold Standard' for diagnosing equine airway pathology. The Endoscope is lightweight and halter mounted and can be quickly moved from one horse to the next. There is nothing that the rider needs to do except complete the normal workout!

Once the insertion tube is placed into the larynx, the distal tip is articulated with the control body. The placement gives an excellent view of the epiglottis and larynx, both, at rest and during the workout. The user can watch the exam 'live' with either wifi or a cellular connection, as well as record the audio and video to the phone. You can then easily transfer the video to a computer with a USB cable, through wifi or with a cellular network. The Tele-View® Dynamic Endoscope comes in two models with different insertion tube diameters and different flexibility.





TELE-VIEW® DYNAMIC  
**EQUINE EXERCISE ENDOSCOPE™**

**Exercise Endoscope Specifications**

**Harness Weight: 4.4 lbs (2.0 KG)**

**Insertion Tube Diameter: 11.5mm or 9.8 mm**

**Endoscope Insertion Tube Length: 42" (106 cm), Custom Lengths Available**

**Resolution: 1280 x 720**

**Recording Camera: Android Cell Phone**

**Cell Phone Charger: 5v USB Phone Charger**

**Working Temperature: 32° - 140° F (0° -60°C)**

**Video Screen: 6.4" (16.2cm) High Resolution**

**Connections: USB C**

**Functions: Real-Time Imaging, Record, Snapshot, Audio**

**Video Output Format: MP4**

*"So yesterday I used your scope the first time on a racing Standardbred. It worked out perfectly! The horse supported very well the scope and the harness. The most positive thing was that the image was very, very stable, even at 1 min 25sec for 1000m the image was perfect. I am very, very glad, that I chose this scope and I would like to congratulate you and your team for this great diagnostic tool."*

*- Dr. Olivier Brandenburg, Dr's. Leser, Korner & Brandenburg, Tierklinik Wiugrund, Rotenburg, Germany*

*"I have used the scope 6 times; all has been smooth sailing and I have been very happy with the image quality."*

*- Dr. Mary-Jane Stutsel, Stutsel Equine, Londonderry, NSW, Australia*

*"Scope well received. Already been used several times. We are very happy with it! Thank you very much."*

*- Cam, Head Equine Veterinary Nurse, West Coast Veterinary Hospital, Perth, Australia*

*"Used it on the first 2 year old race horse to diagnose entrapment today. Worked beautiful. Lots of people on hand, of course, and lots of pressure I might add, including their vet. I couldn't be happier!! Sending the horse to the surgery center today."*

*- Dr. J. W. Clement, Lake Hamilton Equine Associates, Hot Springs, AR*

*"I have the dynamic scope working well, using Anydesk app over the cellular network (to view the live images)"*

*- Dr. Dan Humphrey, Horsemesa Hills Equine Clinic, Echunga, SA, Australia*

**Included Components**



**TV-506B Model**




**Normal Horse Larynx**



**Left Laryngeal Hemiplegia**



**Hard-Sided Case**

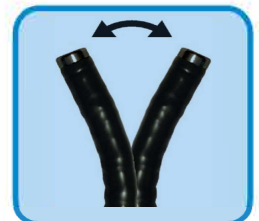

**Articulating Tip**